

# Precision Air Compressor

Air Compressor for Precision Range


Part number: 65-0555 (230V) | 65-1555 (110V)



## Description

The Precision Air Compressor is designed and engineered specifically to provide a dedicated, fully integrated compressed air source for Precision Nitrogen and Precision Zero Air gas generators. Housed within the same compact, stackable form factor as the generators, the Precision Air Compressor module can be added to any Precision stack - typically as the base unit - to give you a complete self-contained GC gas solution, even where no house supply is available.

Thanks to a unique duplex compressor mounting system, the Precision Air Compressor also benefits from low noise and vibration making it relatively quiet in operation within a typical laboratory.

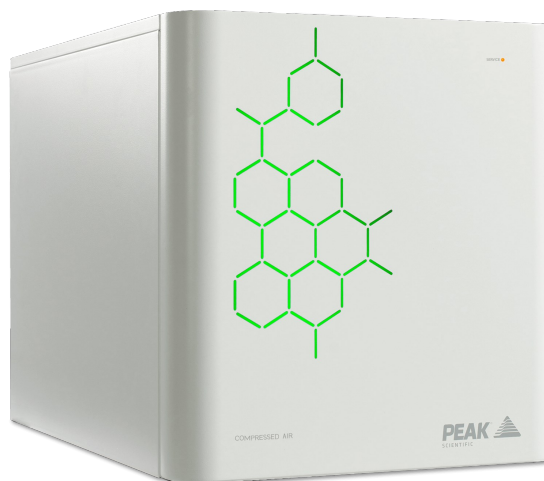

## Application Types


Independent air supply for Precision modules

## Generator Details

- **Gas types:** Air
- **Integrated compressor:** Yes

## Key Features

- Compressed air supply for Precision Nitrogen and Zero Air generators
- Single unit can supply both the Precision Nitrogen and Zero Air generators in a stack
- Duplex compressor mounting system to dampen noise and vibration
- Service & repair status indicators
- 12 month comprehensive on-site warranty

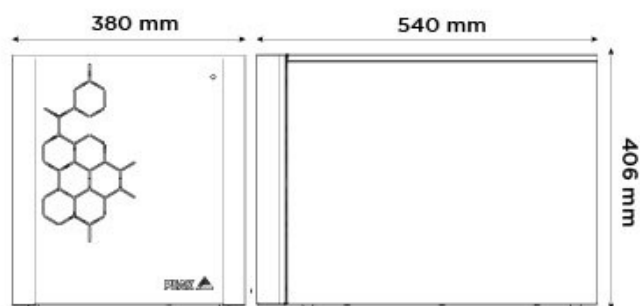

## Technical Specifications

	Precision Air Compressor
Gas Type	Air
Max Gas Flow:	VariousL/min
Max Output Pressure:	120psi/8.3bar
Max Output Pressure:	120psi/8.3bar
Gas Outlets Fitting:	1 x 1/4" BSPP
Power Consumption:	510 / 1150 W
Voltage:	110 / 230 V
Frequency:	50 / 60 Hz
Current:	10 A (110V) / 2.45 A (230V)
Heat output:	1740 - 3922 (BTU)
Max Operating Temp:	35°C / 95°F
Accreditations:	CSA, FCC, CE
Size (HxWxD) mm	406 x 380 x 540 mm
Size (HxWxD) Inches	16 x 15 x 21.3 inches
Generator Weight	42kg / 92.4lbs

## [PEAK Protected]<sup>TM</sup>

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: [www.peakscientific.com/ordering/](http://www.peakscientific.com/ordering/)
- For service plans visit: [www.peakscientific.com/service/service-plans/](http://www.peakscientific.com/service/service-plans/)

**Peak Scientific UK**  
Tel: +44 (0)141 812 8100  
Fax: +44 (0)141 812 8200

**Peak Scientific North America**  
Tel: +1-800-767-6532  
Fax: +1-978-608-9503

**Peak Scientific China**  
Tel: +86 21 5079 1190  
Fax: +86 21 5079 1191  
For a full list of our worldwide office locations, please visit:  
Web: [www.peakscientific.com](http://www.peakscientific.com)  
Email: [discover@peakscientific.com](mailto:discover@peakscientific.com)

Peak Scientific's Quality Management System conforms to: ISO:9001

