

TOC 1500 HP Total Organic Carbon Generator

Application: TOC Analyzers

Part number: 42-1040 (230V) - UK | 42-1050 (230V) - EU | 42-1060 (110V) - US




Description

Peak Scientific's TOC range of specialist generators have been designed specifically to supply clean, dry, carbon dioxide and methane free air for Total Organic Carbon analyzers.

Delivering carrier or combustion gas from an external air supply, these generators are available in two models, depending on your requirement, the [standard model](#) is suitable for non-purgeable organic carbon (NPOC) whilst the HP model is suitable for TC-IC process Total Organic Carbon analyzers.

- **Gas types** Air
- **Integrated compressor:** No
- **Application Types:** TOC Analyzers

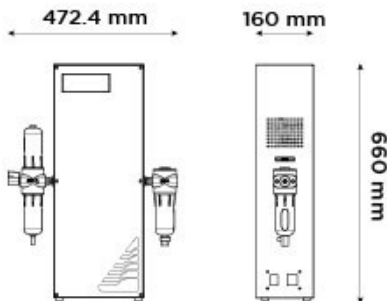

Key Features

- Produces clean, dry, hydrocarbon and carbon dioxide free air
- Ideal for TC-IC process TOC analyzers
- Performs 24/7 if needed
- Few moving parts means little maintenance required and ensures long life of the generator
- Completely silent in operation
- Minimal set-up required
- Can be supported by Precision Air Compressor
- 12 month comprehensive on-site warranty

Technical Specifications

TOC 1500 HP Total Organic Carbon Generator	
Gas Type	Air
Max Gas Flow:	1.5L/min
Max Output Pressure:	80psi/5.5bar
Max Output Pressure:	80psi/5.5bar
Max Purity:	<1ppm CO2/ 0.1ppm THC
Gas Outlets Fitting:	1x 1/4" BSP Female
Min Inlet Pressure:	100psi/6.9bar
Max Inlet Pressure:	120psi/8.3bar
Power Consumption:	115 W
Voltage:	115 / 230 V ± 10%
Frequency:	50 / 60 Hz
Current:	2.6 / 5.2 A
Max Operating Temp:	25°C / 77°F
Accreditations:	CE
Size (HxWxD) mm	620 x 430 x 410 mm
Size (HxWxD) Inches	24.4 x 16.9 x 16.1 inches
Generator Weight	57kg / 125.4lbs

[PEAK Protected]TM

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: www.peakscientific.com/ordering/
- For service plans visit: www.peakscientific.com/service/service-plans/

Peak Scientific UK
Tel: +44 (0)141 812 8100
Fax: +44 (0)141 812 8200

Peak Scientific North America
Tel: +1-800-767-6532
Fax: +1-978-608-9503

Peak Scientific China
Tel: +86 21 5079 1190
Fax: +86 21 5079 1191
For a full list of our worldwide office locations, please visit:
Web: www.peakscientific.com
Email: discover@peakscientific.com

Peak Scientific's Quality Management System conforms to: ISO:9001



0003