

Infinity 1046 Nitrogen Generator

Application: Nitrogen & Air Generator for PerkinElmer QSiight Triple Quad Single Source or Dual Source


Part number: 3302084



Description

Generating up to 20 L/min nitrogen and 70 L/min dry air, the Infinity 1046 has been specifically designed to support the PerkinElmer QSiight Triple Quad, whether single/dual source configuration. The Infinity Series are compressor-free gas generators which require a clean source of in-house air to operate effectively. With very few moving parts, these generators are effectively silent in operation and need only minimal care to ensure continued high performance and extended product life.

- **Gas types** Nitrogen, Air
- **Integrated compressor:** No
- **Application Types:**
Nitrogen & Air Generator for PerkinElmer QSiight Triple Quad Single Source or Dual Source

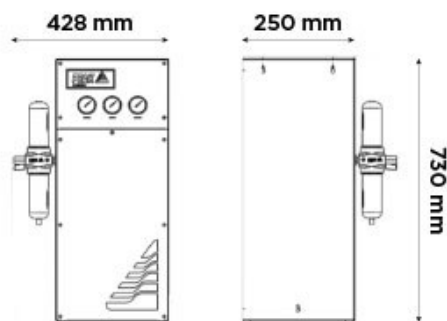

Key Features

- Delivers nitrogen (drying gas) at up to 20 L/min and air (source gas) at up to 70 L/min
- Capable of supplying both single and dual source QSiight models
- 24/7 operation at optimum performance if required
- Highly economical source of nitrogen and air with low lifetime running costs
- Minimal set-up required
- Completely silent in operation
- 12 month comprehensive on-site warranty

Technical Specifications

Infinity 1046 Nitrogen Generator	
Gas Type	Nitrogen
Max Gas Flow:	20L/min
Max Output Pressure:	110psi/7.6bar
Gas Type	Air
Max Gas Flow:	70L/min
Max Output Pressure:	110psi/7.6bar
Max Output Pressure:	110psi/7.6bar
Gas Outlets Fitting:	2 x 1/4" BSPP
Min Inlet Pressure:	120psi/8.3bar
Max Inlet Pressure:	145psi/10bar
Min Air Inlet Flow:	180
Min Inlet Air Quality:	ISO8573 - 1:2010 Class 1.4.1
Start Up Time:	60 mins
Power Consumption:	0
Voltage:	n/a
Frequency:	n/a
Current:	n/a
Heat output:	n/a
Max Operating Temp:	30°C / 86°F
Particles:	<0.01µm
Accreditations:	CE
Size (HxWxD) mm	730 x 422 x 250 mm
Size (HxWxD) Inches	28.7 x 16.6 x 9.8 inches
Generator Weight	23kg / 50.6lbs

.Notes here

[PEAK Protected]TM

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: www.peakscientific.com/ordering/
- For service plans visit: www.peakscientific.com/service/service-plans/

Peak Scientific UK

Tel: +44 (0)141 812 8100
Fax: +44 (0)141 812 8200

Peak Scientific North America

Tel: +1-800-767-6532
Fax: +1-978-608-9503

Peak Scientific China

Tel: +86 21 5079 1190
Fax: +86 21 5079 1191

For a full list of our worldwide office locations, please visit:
Web: www.peakscientific.com
Email: discover@peakscientific.com

Peak Scientific's Quality Management System conforms to: ISO:9001



0003