

NG5000A ultra high purity nitrogen generator

Application: GC, ICP, DSC, DMA

Part number: 10-1292 (230v) - UK | 10-1295 (230v) - EU | 10-1291 (110v) - US




Description

The NG5000(A) (CMS) Carbon Molecular Sieve Generator delivers ultra high purity nitrogen at trace detection limits through Pressure Swing Adsorption technology at flows of up to 5 L/min for GC, DSC, DMA and other applications requiring ultra high purity nitrogen at up to 5L/min. The NG5000(A) comes with an internal compressor for laboratories without an external source of gas.

To reach their ultra high purity nitrogen output these generators use Pressure Swing Adsorption and Carbon Molecular Sieve. These separation technologies remove oxygen and other impurities in air, delivering hassle-free, ultra high purity nitrogen, on-demand for the lab.

- **Gas types** Nitrogen
- **Integrated compressor:** Yes
- **Application Types:** GC, ICP, DSC, DMA

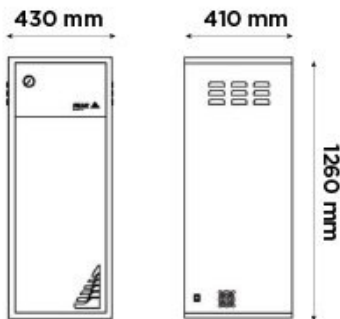

Key Features

- Suitable for applications such as GC, DSC, DMA Suitable for applications requiring ultra high purity
- Available with or without an external compressor
- Contains self-regenerating Carbon Molecular Sieve column
- Gas is supplied on demand so generator works to your schedule
- 12 month comprehensive on-site warranty

Technical Specifications

NG5000A ultra high purity nitrogen generator	
Gas Type	Nitrogen
Max Gas Flow:	5L/min
Max Output Pressure:	80psi/5.5bar
Max Output Pressure:	80psi/5.5bar
Max Purity:	upto 99.9995%
Min Inlet Pressure:	120psi/8.3bar
Max Inlet Pressure:	130psi/9bar
Start Up Time:	8 hours
Power Consumption:	132 Watts
Voltage:	110 / 230 ± 10% VAC
Frequency:	50 / 60 Hz
Current:	0.5 Amps / 1.2Amps
Max Operating Temp:	25°C / 77°F
Accreditations:	CE
Size (HxWxD) mm	1235 x 449 x 406.5 mm
Size (HxWxD) Inches	48.6 x 17.7 x 16 inches
Generator Weight	88kg / 193.6lbs

[PEAK Protected]TM

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: www.peakscientific.com/ordering/
- For service plans visit: www.peakscientific.com/service/service-plans/

Peak Scientific UK
Tel: +44 (0)141 812 8100
Fax: +44 (0)141 812 8200

Peak Scientific North America
Tel: +1-800-767-6532
Fax: +1-978-608-9503

Peak Scientific China
Tel: +86 21 5079 1190
Fax: +86 21 5079 1191
For a full list of our worldwide office locations, please visit:
Web: www.peakscientific.com
Email: discover@peakscientific.com

Peak Scientific's Quality Management System conforms to: ISO:9001



0003