

## NG4000 Nitrogen Generator


Part number: 10-2260 (230v) - UK | 10-2270  
(230v) - EU | 10-2280 (110v) - US



### Description

The NG4000 delivers ultra high purity nitrogen at trace detection limits through Pressure Swing Adsorption technology at flows of 4 L/min for GC applications.

- **Gas types** Nitrogen
- **Integrated compressor:** No

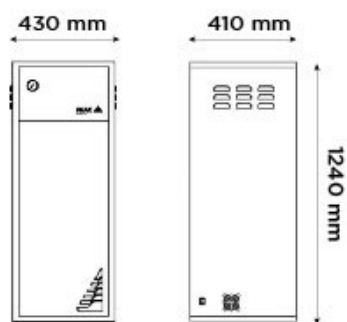

### Key Features

- Suitable for applications such as GC, DSC, DMA
- Suitable for applications requiring ultra high purity
- Available with or without an external compressor
- Contains self-regenerating Carbon Molecular Sieve column
- Gas is supplied on demand so generator works to your schedule
- 12 month comprehensive on-site warranty

## Technical Specifications

	NG4000 Nitrogen Generator
Gas Type	Nitrogen
Max Gas Flow:	40000cc/min
Max Purity:	99.9995%
Min Inlet Pressure:	120 psi
Max Inlet Pressure:	130 psi
Min Air Inlet Flow:	4000 cc/min
Start Up Time:	8 Hours
Voltage:	110/230
Frequency:	50/60
Current:	0.5 Amps
Max Operating Temp:	25°C / 77°F
Particles:	<0.01 µm
Noise Level:	0
Size (HxWxD) mm	1222 x 432 x 406 mm
Size (HxWxD) Inches	48.1 x 17 x 16 inches
Generator Weight	88kg / 193.6lbs

## [PEAK Protected]<sup>TM</sup>

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: [www.peakscientific.com/ordering/](http://www.peakscientific.com/ordering/)
- For service plans visit: [www.peakscientific.com/service/service-plans/](http://www.peakscientific.com/service/service-plans/)

**Peak Scientific UK**  
Tel: +44 (0)141 812 8100  
Fax: +44 (0)141 812 8200

**Peak Scientific North America**  
Tel: +1-800-767-6532  
Fax: +1-978-608-9503

**Peak Scientific China**  
Tel: +86 21 5079 1190  
Fax: +86 21 5079 1191  
For a full list of our worldwide office locations, please visit:  
Web: [www.peakscientific.com](http://www.peakscientific.com)  
Email: [discover@peakscientific.com](mailto:discover@peakscientific.com)

Peak Scientific's Quality Management System conforms to: ISO:9001



0003